



ATLAS MEDICAL

The Health Group

INTRODUCERS, BOUGIES, STYLETS AND GUIDES

We have range of Endotracheal tube Bougies and Stylets, providing clinicians with a comprehensive range of airway devices.

Tracheal Tube Introducer (Bougies) *Straight tip, Curved tip, Vented*

FEATURES AND BENEFIT:

- The Introducer / Bougies used to facilitate endotracheal intubation.
- Made of medical grade PE, LATEX Free, Flexible.
- Curved p ãcan be passed blindly into the trachea when visualizaon ãof the glos ãis inadequate to guide placement of endotracheal tube.
- Distal rounded smooth p ãallows the bougies to be inserted without damaging the so ãssue.ã
- Memory retenon allowãs shaping to individual preferences and needs.
- Clearly graduated marking provides depth indicaon tão ensure correct placement.
- Have different sizes, suitable for pediatric, child and adult.
- Any other special design for your requirements.

PVC Free

DEHP Free

LATEX Free

- Tip types:


Straight tip



Curved tip



Vented tip

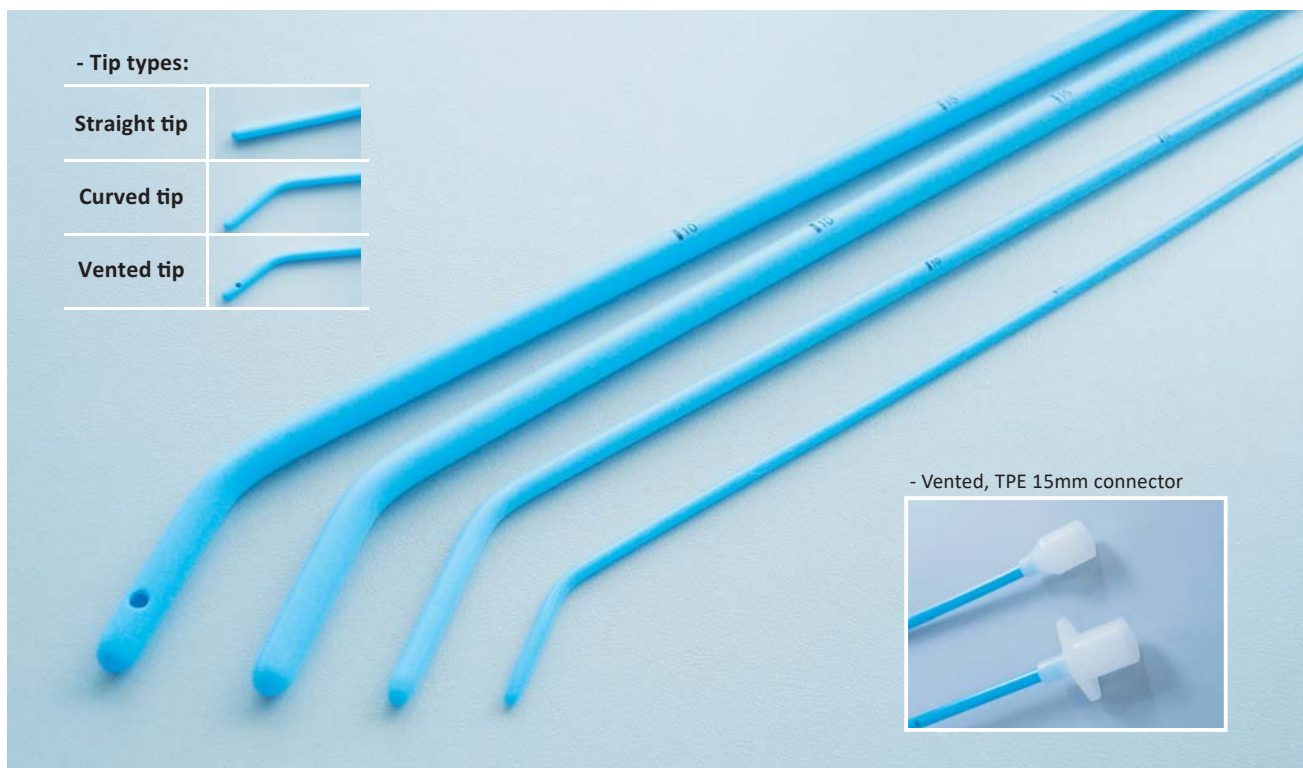

- Vented, TPE 15mm connector

* We can do different color Bougies special for you:



Light Blue



Yellow



Green



Orange



Blue




Purple

INTRODUCERS, BOUGIES, STYLETS AND GUIDES


We have range of Endotracheal tube Bougies and Stylets, providing clinicians with a comprehensive range of airway devices.

1. Tracheal Tube Introducer (Bougies), Curved tip

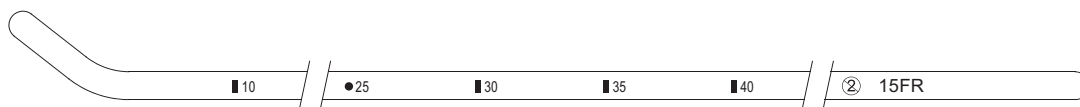


REF number	Size	OD	Length	Recommended Tracheal tube sizes ID (mm)
PA210106	6Fr	2.0mm	53.5cm	2.5 - 3.5
PA210110	10Fr	3.3mm	70cm	4.0 - 5.5
PA210115	15Fr	5.0mm	70cm	6.0 - 10.0

2. Tracheal Tube Introducer (Bougies), Straight tip

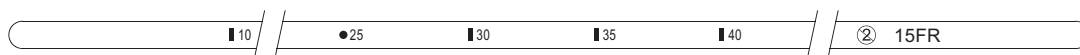

REF number	Size	OD	Length	Recommended Tracheal tube sizes ID (mm)
PA210206	6Fr	2.0mm	53.5cm	2.5 - 3.5
PA210210	10Fr	3.3mm	70cm	4.0 - 5.5
PA210215	15Fr	5.0mm	70cm	6.0 - 10.0

3. Tracheal Tube Introducer (Bougies), Vented tip

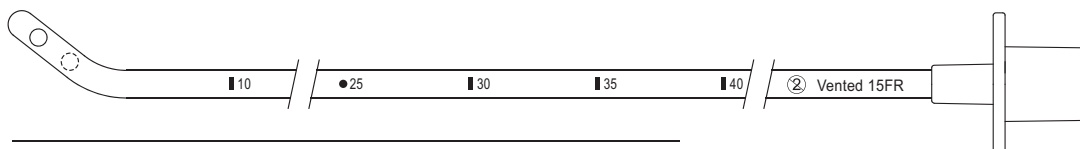

REF number	Size	OD	Length	Recommended Tracheal tube sizes ID (mm)
PA210308	8Fr	2.7mm	53.5cm	3.0 - 4.0
PA210310	10Fr	3.3mm	70cm	4.0 - 5.5
PA210314	14Fr	4.7mm	70cm	6.0 - 10.0
PA210315	15Fr	5.0mm	70cm	6.0 - 10.0

4. Tracheal Tube Introducer (Bougies), Curved tip, Easy-Carry Pack

FEATURES: InnovavÆ compact packaging designed specifically for use by EMS professionals.

REF number	Size	OD	Length	Recommended Tracheal tube sizes ID (mm)
PA210510	10Fr	3.3mm	70cm	4.0 - 5.5
PA210515	15Fr	5.0mm	70cm	6.0 - 10.0

