

SINGLE USE LINEAR CUTTING
STAPLER AND RELOAD



ATLAS MEDICAL

The Health Group

OPEN INSTRUMENTS

TCYQ Series




TONCARE
Medical Instrument
Manufacture Center

CE 0197 / ISO 13485


ADVANCING FOR LIFE

SINGLE USE LINEAR CUTTING STAPLER AND RELOAD

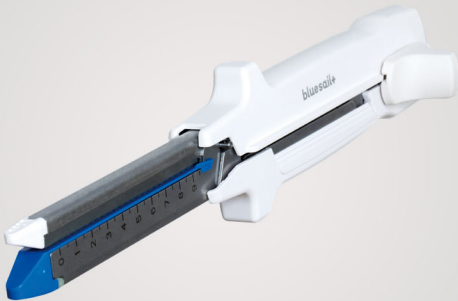

Product Overview

TCYH Series has application in gastrointestinal, gynecologic, thoracic and pediatric surgery for transection, resection and/or creation of anastomoses.

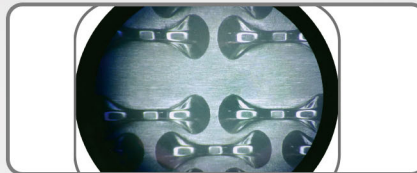

Droplet Groove Design

Increase oriented contact area to ensure better B-shape staple formation



CAM Design

Provide consistent formation from proximal to distal to reduce risk of staple line leakage



Safety Lock-out Design

Prevent a used reload from being refired



High Strength Staple

Use High strength titanium staple to get better formation

Product Code	Color	Stapling Length (mm)	Number of Staples	Tissue Thickness (mm)	Staple Height (mm)
TCYQ-60×3.85	Blue	60	60	1.5	3.85
TCYQ-60×4.5	Green	60	60	2.0	4.50
TCYQ-75×3.85	Blue	75	75	1.5	3.85
TCYQ-75×4.5	Green	75	75	2.0	4.50
TCYQ-90×3.85	Blue	90	90	1.5	3.85
TCYQ-90×4.5	Green	90	90	2.0	4.50
TCYQ-100×3.85	Blue	100	100	1.5	3.85
TCYQ-100×4.5	Green	100	100	2.0	4.50



ATLAS MEDICAL

The Health Group



Gallagher Place North U-14
609 Richards Road
Halfway House X94
Midrand
1685



+27 (11) 312 0096



www.atlasmedical.co.za



sales@atlasmedical.co.za

SINGLE USE LINEAR
CUTTING STAPLER
AND RELOAD